

Car Extractor

The Car Extractor is designed to safely lift and remove a patient from nearly any vehicle, from a compact automobile to a fire/emergency van. It is ideal for hospital emergency departments as this heavy-duty lift enables speedy transfers of patients onto gurneys, stretchers or wheelchairs. The car extractor can help insert patients into vehicles effectively, too.


Built light and strong, with a 30" reach and a 40" lifting range, the car extractor can gain deep access into tight vehicle interior spaces and lift patients with limited mobility or little weight bearing ability—up to 475 lbs. capacity.

Item Number:

450000 – Car Extractor



“Perfect for reducing injuries related to vehicle transfers, the Car Extractor is ideal for transferring patients into and out of most vehicles.”

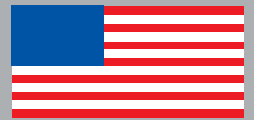

MADE IN THE USA



Car Extractor

Item #450000 – lifts up to 475 lbs.



Technical Specifications

Overall Length:	47"
Outside Width (legs closed):	24"
Inside Width (legs open)	46"
Height of Wheel Base	4.5"

Boom Arm Reach	30"
Lifting Range	40"
Power Source	Rechargeable 24 volt, 5 amp battery

Med
HANDLE WITH CARE.

Accessories:

Hanger Bar Stopper Option

The Medicare Hanger Bar Stopper prevents sling loops from sliding off of the hanger bar hooks.



Item Number:

400611

PRODUCT HIGHLIGHTS

- Dual battery design allows lift to be in constant use
- 360 degree rotation for easy positioning
- Straight boom for deep access into vehicle
- Electronic actuator lifting system — no hydraulic lifting
- Remote hand control to allow operator to maintain contact with patient
- Modular design ensures lift will never be obsolete and can always be upgraded
- Auxiliary switch for raising and lowering lift
- Compact, padded 2-point hanger bar
- Low leg base will slide under almost any vehicle
- To be used with two Medicare Stand Belts and two personnel



Medcare Products

151 East Cliff Road
Burnsville, MN 55337

Toll-free: 800-695-4479

Fax: 952-894-7153

www.MedicareLifts.com

© 2011 Medicare Products Inc.
Specifications are subject to change without notice. All products referred herein are registered trademarks or trademarks of their respective holders.