

Medcare Care Lifts are designed to be used by caregivers to transfer non-weight bearing patients without the strain of manual transfers. Available in two weight capacities, 475 lbs. and 600 lbs., which is our bariatric capacity lift. Medcare Lifts eliminate much of the exposure to transfer-related injuries. The Care Lift-N-Weigh[®] and Care Bâr Lift-N-Weigh[®] feature our convenient scale option. The expandable base allows for easy access around wheelchairs, beds, tubs and toilets; creating a safe and user-friendly transfer system.

Item Number:

400001 – Medcare Care Lift, lifts up to 475 lbs.

400019 – Care Bâr Lift, lifts up to 600 lbs.

400003 – Care Lift-N-Weigh (includes scale), lifts up to 475 lbs.

400013 – Care Bâr Lift-N-Weigh (includes scale), lifts up to 600 lbs.




Care Lift



Lift-N-Weigh

“Lifts are an important part of Medcare’s zbi[®] (Zero Back Injury) system which reduces costs associated with workers-compensation claims and medical expenses.”

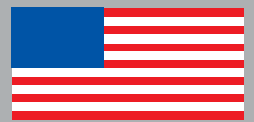

MADE IN THE USA

Features:

- Easily transfers patients to and from chairs, wheelchairs, and toilets
- Floor and bed transfers can be accomplished with ease
- Designed to work with a variety of Medcare slings to “customize” the transfer





Hand Control
Item # 400990



Battery Charger

Battery

Med
HANDLE WITH CARE.

Accessories:

Hanger Bar Stopper Option

The Medicare Hanger Bar Stopper prevents sling loops from sliding off of the hanger bar hooks.



Item

Number:

400611

Scale Option

Medcare's scale option, available on the Care Lift-N-Weigh and Care Bar Lift-N-Weigh, provides a quick efficient means of transferring and weighing patients in one easy step. The scale's display is mounted on the mast so that caregivers can easily read and record a patient's weight.



Technical Specifications for 475 lbs. & 600 lbs. Lifts

Overall Length:	47"	Height of Wheel Base	4.5"
Outside Width (legs closed):	24"	Power Source	Rechargeable 24 volt, 5 amp battery
Inside Width (legs open)	46"		

PRODUCT HIGHLIGHTS

- Dual battery design allows lift to be in constant use
- 360° rotation for easy patient positioning
- Arched boom for added head clearance
- Electronic actuator lifting system — no hydraulic lifting
- 4-point, non-rigid hanger bar design allows sling to conform to patient for added comfort
- Remote hand control to allow operator to maintain contact with patient
- Modular design ensures lift will never be obsolete and can always be upgraded
- Auxiliary switch for raising and lowering lift
- Low leg base allows use under lower bed rails; 4.5" with standard casters, 3" low-profile casters available
- To be used with Medicare Slings



Medcare Products

151 East Cliff Road
Burnsville, MN 55337

Toll-free: 800-695-4479

Fax: 952-894-7153

www.medcarelifts.com

© 2010 Medcare Products Inc.
Specifications are subject to change without notice. All products referred herein are registered trademarks or trademarks of their respective holders.