

Care Stands

Medcare Care Stands are designed to assist caregivers in raising weight-bearing patients to a standing position without the need for manual lifting. The compact design provides easy access to confined spaces such as bathrooms. Available in two weight capacities, 350 lbs. and a bariatric unit lifting up to 500 lbs.

The Care Stand-N-Weigh® and Care Bâr Stand-N-Weigh® both feature our patented scale.

Item Number:

400002 – Medcare Care Stand, lifts up to 350 lbs.

400006 – Care Bâr Stand, lifts up to 500 lbs.

400005 – Care Stand-N-Weigh (includes scale), lifts up to 350 lbs.

400007 – Care Bâr Stand-N-Weigh (includes scale), lifts up to 500 lbs.




Care Stand



Care Stand-N-Weigh®

“Stands are an important part of Medcare’s **zbi**® (Zero Back Injury) system which reduces costs associated with workers-compensation claims and medical expenses.”

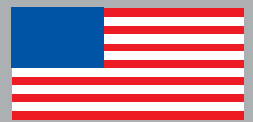

MADE IN THE USA

Features:

- Removable foot plate provides ambulation capability
- Easy toileting
- Enhanced patient hygiene
- Convenient for changing patient’s garments





Hand Control
Item # 400990



Battery Charger

Battery



Technical Specifications for 350 lbs. & 500 lbs. Stands

Overall Length:	39"	Height of Wheel Base	4.5"
Outside Width (legs closed):	24"	Power Source	Rechargeable 24 volt, 5 amp battery
Inside Width (legs open)	43"		

PRODUCT HIGHLIGHTS

- Dual battery design allows stand to be in constant use
- Hand control allows caregiver to be positioned next to patient during transfer
- Padded knee support is adjustable and removable and includes a Velcro® strap to keep patient legs secured if necessary
- Removable floor platform is flat and parallel to the ground so patient can stand upright
- 4.5" leg base height allows use with low beds
- Electronic actuator lifting system — no hydraulic lifting
- Modular design ensures stand will never be obsolete and can always be upgraded
- Auxiliary switch for raising and lowering stand
- To be used with Medcare Belts and Transport Belts
- No need for a separate scale — weigh patients as they are transferred
- Scale can be used with ambulatory patients, just have them step on the foot plate and record weight

Med
HANDLE WITH CARE.

Accessories:

Hanger Bar Stopper Option

The Medcare Hanger Bar Stopper prevents sling loops from sliding off of the hanger bar hooks.



Item

Number:

400611

Scale Option

Medcare's patented scale option, available on the Care Stand-N-Weigh® and Care Bâr Stand-N-Weigh®, allows standing transfers to be



performed while recording a patient's weight. Whether a patient is standing straight, at an angle or even holding onto the handles, the scale accurately weighs patients.



Medcare Products

151 East Cliff Road
Burnsville, MN 55337

Toll-free: 800-695-4479

Fax: 952-894-7153

www.MedcareLifts.com

© 2011 Medcare Products Inc.
Specifications are subject to change without notice. All products referred herein are registered trademarks or trademarks of their respective holders.